



SuperCath[®] 5

26G Safety IV Catheter



Safety peripheral IV catheter technology designed to help you access small, difficult vessels.


icumedical
human connections

Helping You Deliver Safe and Effective Care While Complying with Clinical Recommendations

Designed to optimise vessel entry in your neonatal patients, the SuperCath 5 is the first and only 26G straight safety catheter. This unique needle size and integrated needle notch helps you improve first-stick proficiency and avoid multiple insertion attempts by maximising the speed and likelihood of flash visualisation while accessing small, difficult vessels.

Infusion Therapy Standards of Practice recommendation:

Always select the smallest gauge peripheral catheter that will accommodate the prescribed therapy and patient need.¹

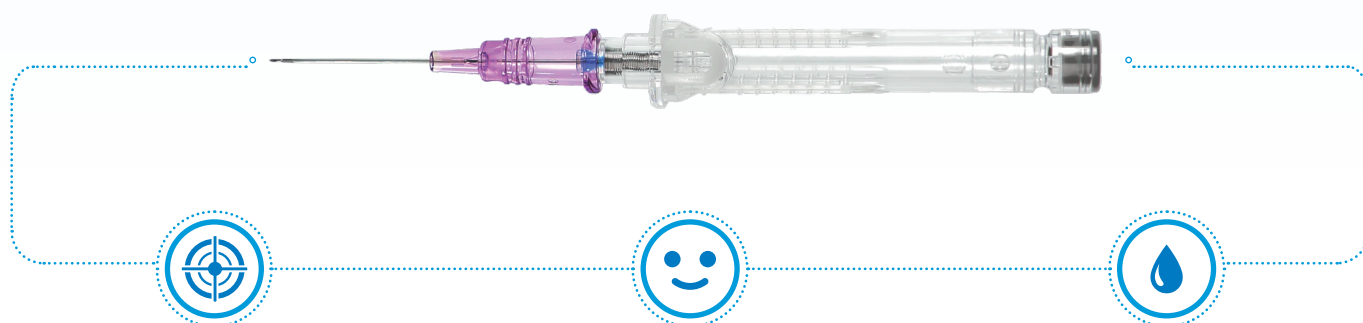

Maximise First Stick Proficiency

Help ensure the catheter is properly in the vessel on the first stick with needle notch technology.

Minimise Patient Discomfort

Help reduce patient discomfort and the need for multiple insertion attempts with small (26G) needle with notch.

Reduce the Risk of Needlestick Injury

Avoid accidental needlesticks and needle exposure with push button, spring-loaded needle encapsulation.

Product Information

Product Code Non-Winged	Product Code Winged	Gauge	Length		O.D. (mm)	I.D. (mm)	Flow Rate (mL/min)
			(in)	(mm)			
SP125-26-14T*	SP125-26-14-WT*	● 26	0.56	14	0.6	0.4	14
SP125-26-19T*		● 26	0.75	19	0.6	0.4	12

*No blood-containment valve.

1. Infusion Nurses Society. Infusion nursing standards of practice. J Infus Nurs. 2016; 39 (suppl 1): S51

SuperCath is a registered trademark of MediKit Co., Ltd.

The product complies with current legislation and has the corresponding CE marking. For additional information, warnings and / or safety precautions, refer to the manufacturer's user manual.

© 2023 ICU Medical, Inc. | P23-4630-EMEA

icumedical
human connections

+44 (0)203 357 9400 | www.icumed.com