



Plum 360™ Infusion System Catalogue



Connect Smarter. Care Safer™

Get award-winning clinical performance.

The Best in KLAS Plum 360 infusion system gives you automatic safety with unmatched efficiency and interoperability in one connected IV system.

The Plum 360 features a unique cassette and delivery system which enables

- › A direct connection of the secondary line to the Plum™ cassette, so the pump knows which container—primary or secondary—it is delivering from
- › Concurrent delivery of two compatible medications at independent rates through a single IV line
- › Capture of 1 mL of air before alarming, helping to reduce air-in-line alarms
- › Integration with your EHR for enhanced safety and efficiency



Order Code	Item Number
443001007	300105421
Description	
Plum 360 Infuser compatible with ICU Medical MedNet™ software	




Best in KLAS
8 years
in a row!



Dedicated Plum Sets

The Plum cassette is at the heart of the unique Plum delivery system. The direct connection of the secondary line to the cassette provides control and management for both primary and secondary infusions. You can feel confident that both your piggyback and concurrent infusions (of compatible medications) are delivering as programmed.

The Plum cassette includes an air trap that captures up to 1 mL of air before alarming. This has been shown to reduce air-in-line alarms. If an air-in-line alarm does occur, air can be easily removed without disconnecting from the patient. This closed system also helps to reduce the risk of contamination and limit exposure to hazardous medications.

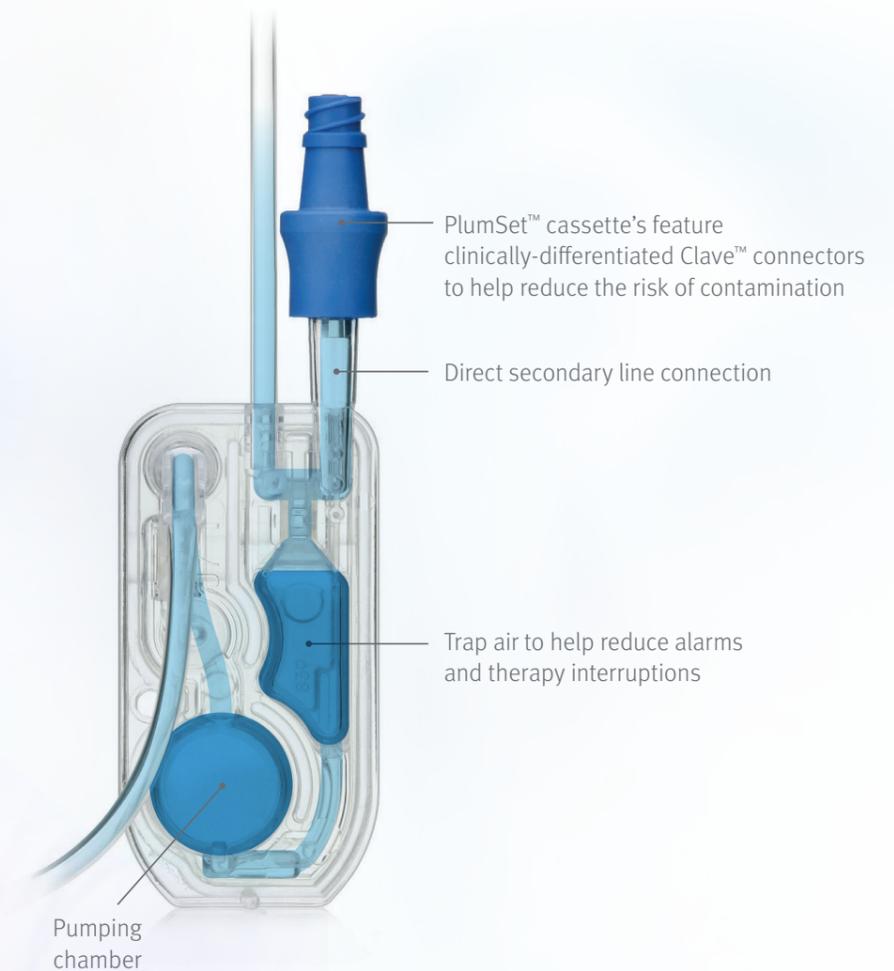


Plum 360 Primary Sets


Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441400004	14000		Primary Plum Set, 15 Micron Filter in Sight Chamber, Prepierced Y-Site, Secure Lock, 272 cm /19 mL	272	19.0	PVC non-DEHP	50




Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441400104	14001	FSB1314	Primary Plum Set, 15 Micron Filter in Sight Chamber, Clave Secondary Port, Clave Y-Site, Secure Lock, 272 cm/19 mL	272	19.0	PVC non-DEHP	50

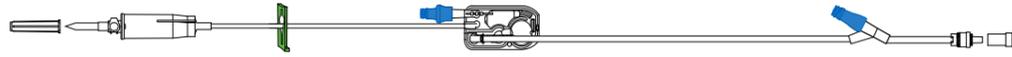


Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
	14003		Primary Plum 150 mL Burette Set, Prepierced Additive Port, 15 Micron Filter in Sight Chamber, Prepierced Secondary Port, 2 Prepierced Y-Sites, Secure Lock, 290 cm / 19.4 mL	290	19.4	PVC non DEHP	50

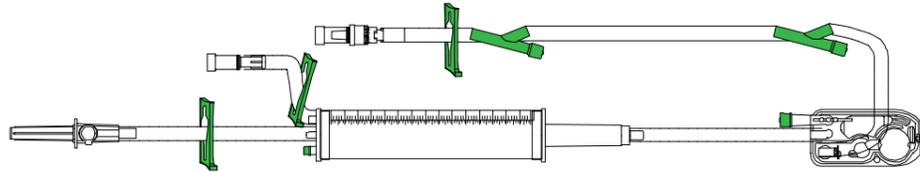

Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441400604	14006		Primary Plum Set, 15 Micron Filter in Sight Chamber, Polyethylene Lined Light Resistant Tubing, Distal Microbore Tubing, Secure Lock, 272 cm/12 mL	272	12.0	PVC non-DEHP/PE	50

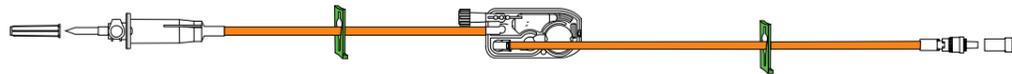


Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441400704	14007	FSB2379	Primary Plum Set, 15 Micron Filter in Sight Chamber, Clave Secondary Port, Clave Y-Site, Polyethylene Lined Light Resistant Tubing, Distal Microbore Tubing, Secure Lock, 272 cm/12 mL	272	12.0	PVC non-DEHP/PE	50

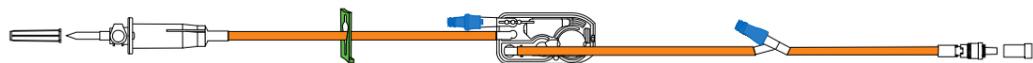


Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441400904	14009	FSB2467	Primary Plum Set, 15 Micron Filter in Sight Chamber, Clave Secondary Port, 0.2 Micron Filter, Clave Y-Site, Polyethylene Lined Light Resistant Tubing, Distal Microbore Tubing, Secure Lock, 272 cm/15 mL	272	15.0	PVC non-DEHP/PE	50




Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441401203	14012	FSB1792	Primary Plum Set, 15 Micron Filter in Sight Chamber, Clave Secondary Port, Secure Lock, 225 cm/17 mL	225	17.0	PVC non-DEHP	50

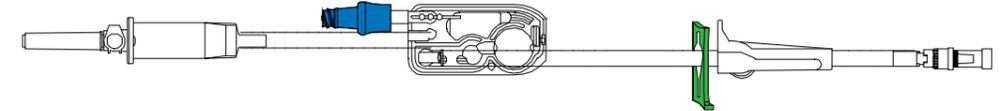

Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441401403	14014	FSB1696	Primary Plum Set, 15 Micron Filter in Sight Chamber, Clave Secondary Port, 0.2 Micron Filter, Polyethylene Lined Tubing, Secure Lock, 225 cm/20 mL	225	20.0	PVC non-DEHP/PE	50

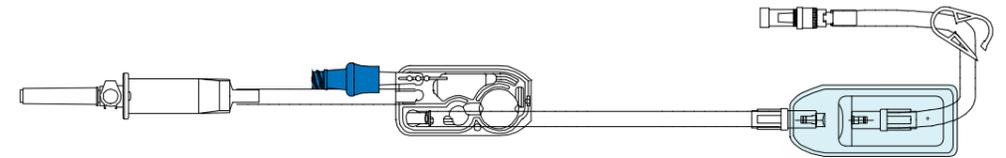


Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441401504	14015		Primary Plum Set, 15 Micron Filter in Sight Chamber, Clave Secondary Port, 0.2 Micron Filter, Clave Y-Site, Polyethylene Lined Tubing, Secure Lock, 272 cm/23 mL	272	23.0	PVC non-DEHP/PE	50

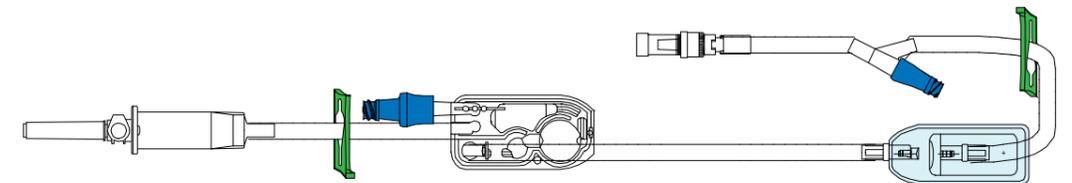


Plum 360 Secondary Sets

With the Plum 360 infusion system's unique primary cassette pumping technology and cassette, you can deliver secondary infusions through a single line with confidence. Setting up a treatment regime with your fluids on the primary line and medications on the secondary line can save valuable time for your patients, and reduce the consumption of flush bags for your facility.

Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
140179290	14017		Primary Plum MICRODRIP Set, 15 Micron Filter in Sight Chamber, Clave Secondary Port, Clave Y-Site, Polyethylene Lined Tubing, Secure Lock, 272 cm/20 mL	272	20.0	PVC non-DEHP/PE	50

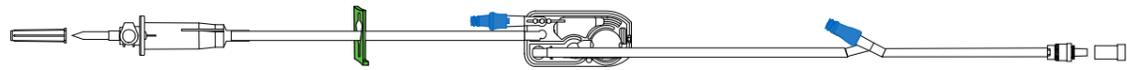


Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441967302	19673	FSB2100	Primary Plum Blood Set, 200 Micron Filter, Clave Secondary Port, Non-Vented, Secure Lock, 259 cm/51 mL	259	51.0	PVC non-DEHP	48

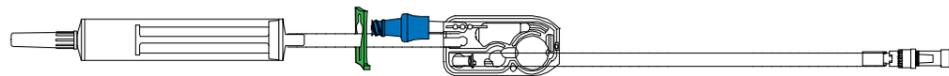

Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441421204	14212	FSB2252	Primary Plum Y-Type Blood Set, 200 Micron Filter, Clave Secondary Port, Non-Vented, Secure Lock, 279 cm/49 mL	279	49.0	PVC non-DEHP	48




Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441256401	12564		Primary Plum 150 mL Burette Set, Clave Additive Port, 15 Micron Filter in Sight Chamber, 2 Clave Y-Sites, Secure Lock, 290 cm/19 mL Retrograde Volume is Approximately: 8.3 mL	290	19.0	PVC non-DEHP	48

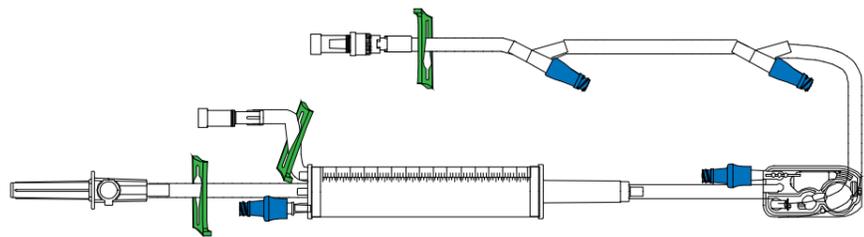


Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
441402804	14028	FSB1398	Secondary Set, 15 Micron Filter in Sight Chamber, Secure Lock, 86 cm / 6 ml	86	6.0	PVC non-DEHP	50



Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
44H326101	011-CH3261	FSB2319	76 cm (30") Appx 3.3 mL, Admin Set, 2 Spiros™ w/Red Cap, 20 Drop In-Line Drip Chamber w/15 Micron Filter, Bag Hanger	76	3.3	PVC non-DEHP	50

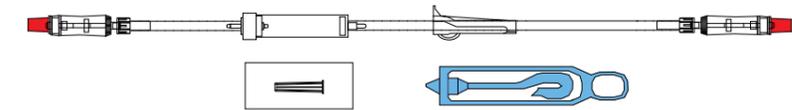


Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
44H323501	CH3235	FSB1240	30" (76 cm) Appx 4.1 mL, Yellow 20 Drop Admin Set w/15 Micron Filter, Spiros™	76	4.1	PUR	50

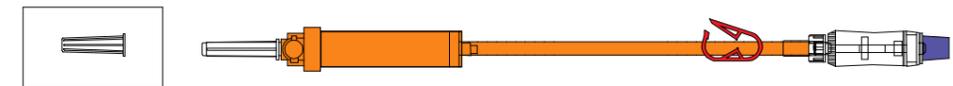

Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
44H330401	011-CH3304	FSB1973	76 cm (30") 20 Drop Admin Set w/Integrated ChemoClave™ Drip Chamber, Bag Hanger, Rotating Luer	76	3.3	PVC non-DEHP	50



Order Code	Item No	NHSSC NPC Code	Description	Length (cm)	PV (mL)	Tubing Material	Case Qty
011-CH3751	011-CH3751	FSB2490	Admin Set, Amber, Spiros	76	6.9	PVC non-DEHP	50

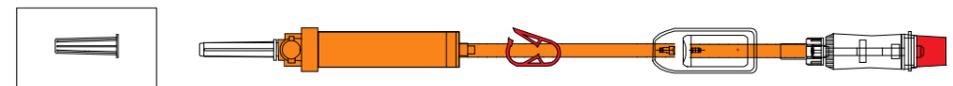

Pump Accessories

ICU Medical provides a range of IV poles and mounting accessories for our Infusion pumps to help you to efficiently organise your space around the patient.



Order Code	Item No	Case Qty
391192502	323840	1
Description		
Provita heavy duty pole. Pole can accommodate up to 5 x Plum 360 pumps. A Triple pump holder and a dual pump holder are required for this configuration		



Order Code	Item No	Case Qty
391273602	323761	1
Description		
Provita premium IV pole. Pole can accommodate up to 4 x Plum 360 pumps. 2 x Dual pump holders or a triple pump holder with a single pump below on the pole		



Order Code	Item No	Case Qty
391052501	175191	1
Description		
Provita standard pole. Pole can accommodate up to 2 x Plum 360 pumps (attached directly to the pole)		

Provita IV Pole Mount Accessories



Order Code	Item No	Case Qty
391233601	324751	1
Description		
Provita dual pump holder. Use to attach two Plum 360 pumps to a premium or heavy duty IV pole		



Order Code	Item No	Case Qty
391202501	324283	1
Description		
Provita triple pump holder. Use to attach three Plum 360 pumps to a premium or heavy duty IV pole		



Order Code	Item No	Case Qty
391303601	324317	1
Description		
Provita 6 point UK plug board—Pin type G. Power multiple pumps with one plug to the wall (3 pin Type G). Universal clamp attaches to poles with different diameters (including Heavy Duty and Premium poles)		



Order Code	Item No	Case Qty
391213602	324292	1
Description		
Provita bag holder. Keep bags handy and lines organised when multiple pumps are mounted to a pole		



Order Code	Item No	Case Qty
391223601	329394	25
Description		
Cable Tidy. Keep cables tidy when multiple pumps are mounted to a pole. Can only be used with Heavy duty pole (391192502)		



Order Code	Item No	Case Qty
391293601	329437	25
Description		
Provita single cable tidy for 25 mm tube Keep your cable tidy on the pole (single cable)		



Order Code	Item No	Case Qty
391283601	329491	1
Description		
Provita pump holder for wall rail 25 x 10 mm Add pumps to a medical rail on the wall in a hospital room		

The product complies with current legislation and has the corresponding CE marking (2797).
For additional information, warnings and/or safety precautions, refer to the manufacturer's Instructions for Use.

icumedical
human connections