



EchoGlo™ Peripheral Nerve Block Solutions

Engineered for enhanced visualization

icumedical
human connections



Engineered for ultrasound




The Portex™ EchoGlo peripheral nerve block needle is designed with two engraved bands that increase ultrasound visibility by optimally reflecting ultrasound waves.



EchoGlo catheters are echogenic from tip to tip

- › The Portex EchoGlo catheter is engineered with flexible, echogenic, nylon tubing.
- › The marked catheter tip provides for ease of identification and to ensure complete catheter removal.

The anatomy of an echogenic needle

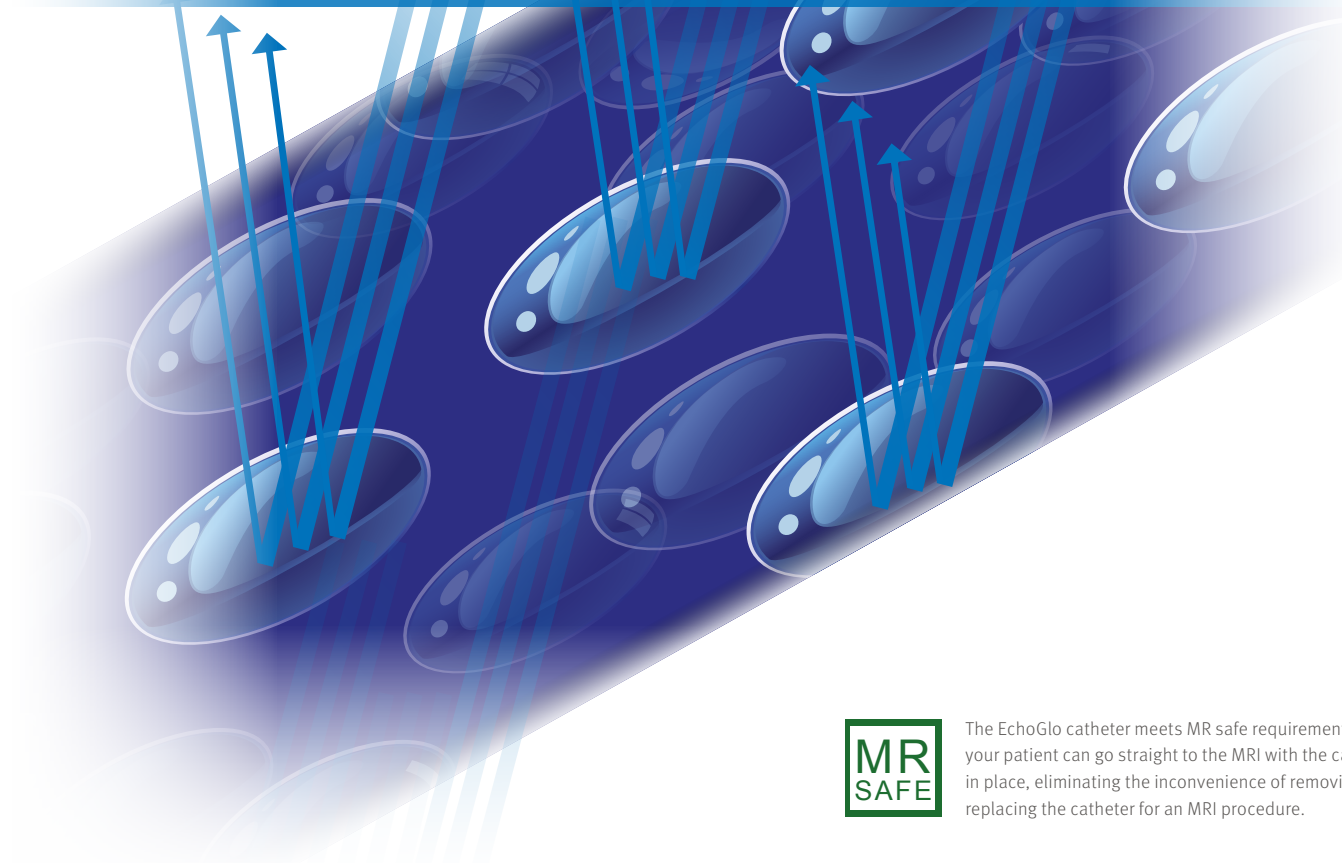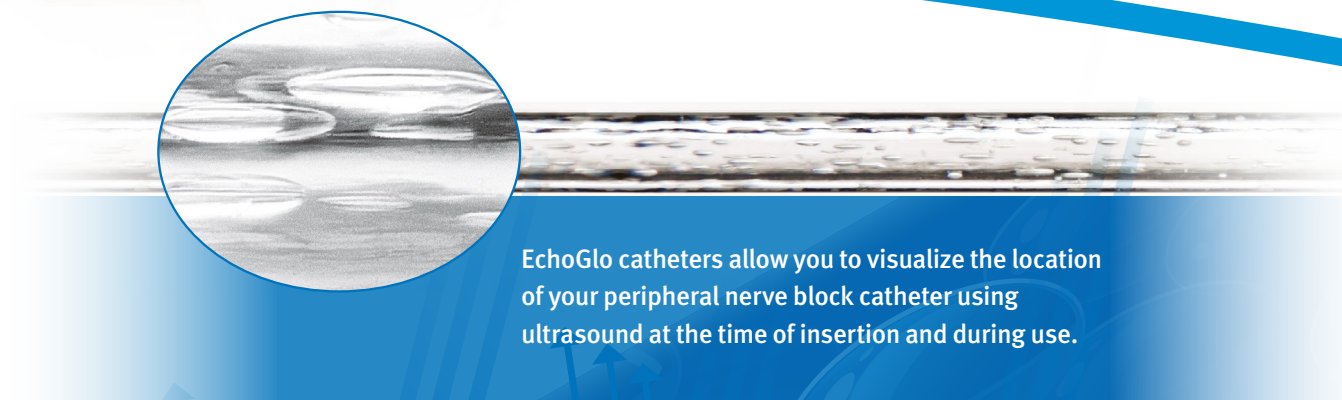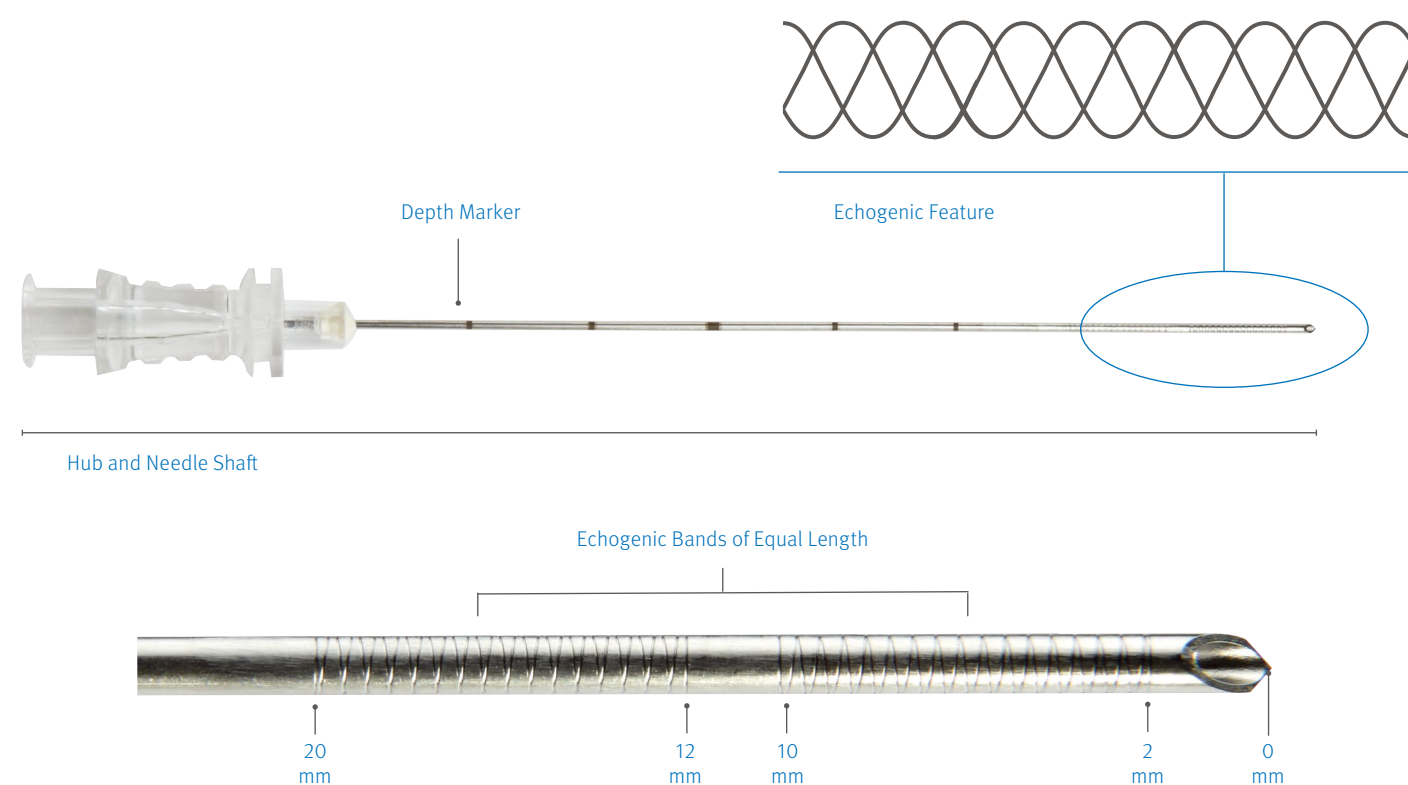

The EchoGlo catheter meets MR safe requirements so your patient can go straight to the MRI with the catheter in place, eliminating the inconvenience of removing and replacing the catheter for an MRI procedure.


Ordering Information

ICU Medical has a full portfolio of Echogenic peripheral nerve block solutions to support both single-shot and continuous procedures.



EchoGlo Peripheral Nerve Block Needle – Single Shot			
Product Number	Description	Metric/Nonmetric	Per Box
24-4100-24	Single-shot nerve block needle, 22G x 35 mm, bevel, Transparent Female Luer	35 mm/1 ¾ in	25
24-4101-24	Single-shot nerve block needle, 22G x 50 mm, bevel, Transparent Female Luer	50 mm/2 in	25
24-4102-24	Single-shot nerve block needle, 22G x 80 mm, bevel, Transparent Female Luer	80 mm/3 ⅛ in	25
24-4103-24	Single-shot nerve block needle, 22G x 100 mm, bevel, Transparent Female Luer	100 mm/4 in	25
24-4104-24	Single-shot nerve block needle, 22G x 100 mm, bevel, Transparent Female Luer	100 mm/4 in	25
24-4105-24	Single-shot nerve block needle, 22G x 150 mm, bevel, Transparent Female Luer	150 mm/6 in	25

Includes a 610 mm (24 in) extension set

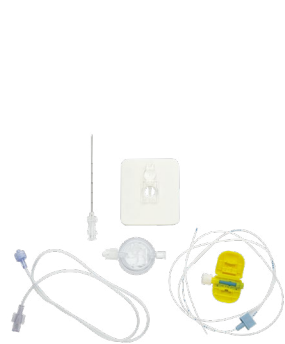

EchoGlo Peripheral Nerve Block Kits – Continuous			
Product Number	Needle Gauge	Metric/Nonmetric	Per Box
24-4114-24	Continuous nerve block kit, 18G x 50 mm, Tuohy	50 mm/2 in	10
24-4115-24	Continuous nerve block kit, 18G x 80 mm, Tuohy	80 mm/3 ⅛ in	10
24-4116-24	Continuous nerve block kit, 18G x 100 mm, Tuohy	100 mm/4 in	10
24-4117-24	Continuous nerve block kit, 18G x 150 mm, Tuohy	150 mm/6 in	10
24-4118-24	Continuous nerve block kit, 18G x 50 mm, Bevel	50 mm/2 in	10
24-4119-24	Continuous nerve block kit, 18G x 80 mm, Bevel	80 mm/3 ⅛ in	10
24-4120-24	Continuous nerve block kit, 18G x 100 mm, Bevel	100 mm/4 in	10
24-4121-24	Continuous nerve block kit, 18G x 150 mm, Bevel	150 mm/6 in	10

Includes continuous needle, 21G closed-end 3 eyest EchoGlo catheter, LockIt Plus™ device, EpiFuse™ catheter connector, extension set, thread-assist guide, and filter.



EchoGlo Peripheral Nerve Block Continuous Trays				
Procedural Needles	24-4126-01	24-4127-01	24-4128-01	Per Case
18G 50 mm (2”) Tuohy	•			10
18G 80 mm (3 1/8”) Tuohy			•	10
18G 100 mm (4”) Tuohy		•		10
Catheter System				
Closed-end EchoGlo catheter	•	•	•	10
EpiFuse™ catheter connector	•	•	•	10
Catheter insertion device	•	•	•	10
LockIt Plus™ catheter securement device	•	•	•	10
Hypodermic Needles				
18G x 1 ½ in	•	•	•	10
22G x 1 ½ in	•	•	•	10
25G x 1 ½ in	•	•	•	10
Syringes				
3 mL Plastic Luer Lock	•	•	•	10
5 mL Plastic Luer Lock	•	•	•	10
20 mL Plastic Luer Lock	•	•	•	10
Other Procedural Components				
Ultrasound probe cover kit	•	•	•	10
24 in extension set	•	•	•	10
Flat Filter	•	•	•	10
Filter Straw	•	•	•	10
Procedural Drugs				
Lidocaine HCl (1%) 5 mL	•	•	•	10
Lidocaine HCl (1.5%) and Epinephrine (1:200,000) 5 mL	•	•	•	10
Sodium Chloride (0.9%) 10 mL	•	•	•	10
Filter Straw	•	•	•	10
Drape				
Clear Fenestrated Drape	•	•	•	10
Other Procedural Components				
Towel	•	•	•	10
ChloraPrep™ One-Step	•	•	•	10
Four Gauze Sponges	•	•	•	10
Needle Stick Pad	•	•	•	10
Nerve Block Catheter I.D. Label	•	•	•	10

